

Standardization of Proton Induced X-Ray Emission for Analysis of Trace Elements

Johar Zeb¹, Shad Ali¹, Abdul Ahad³, Ishfaq Ahmad³, Jehan Akbar²

¹Department of Physics Beijing Normal University, Haidian, China, ²Department of Physics Hazara University Mansehra, ³National center of Physics (NCP), Islamabad (Pakistan)

E-mail: Joharzeb789@yahoo.com

This paper presents the standardization of Proton Induced X-rays Emission (PIXE) technique for the trace element analysis of three standard samples. Standard Reference Materials (SRM's) (Titanium based alloy, Copper based alloy and Iron based alloy) are proposed for the desired PIXE analysis in such a way that it cover the maximum range of elements in the periodic table. Protons beam with three different values charge was used to incident on the sample up to 2.57MeV energy from 5UDH-II tandem Pelletron accelerator for the X-rays emission from the sample. The beam is used to incident at different positions on target materials, each sample is exposed to beam for more than three times for purpose to get X-rays from approximately all the composite elements of the sample. The acquired data was analyzed by GUPIXWIN software for the determination of trace elements. The acquired data machine is compared with SRM data. The comparison between these data was in closed agreement. Finally some discussion is made on our analysis their results in compared with other elemental analytical techniques

Key Words: PIXE Standardization, Trace Elements, GUPIXWIN Software.

1. INTRODUCTION

The word PIXE stands for Particle/Proton Induced X-rays Emission. PIXE was first experimentally shown by Sven A.E. Johansson of Lund Institute of Science and Technology in 1970 [2]. The advantage of this technique is its non-destructive nature and multi-elemental character of its analysis simultaneously. Due to its longer probe depth (tens of μm) the PIXE can be performed for a large range of range of elements ($\text{Na} \leq Z \leq \text{U}$) in the periodic table. PIXE is fast, highly sensitive (in limit ppm), precise, accurate and is better for analysis of both matrix and trace elements in thick sample as compared to other ion beam analytical techniques [2, 3]. PIXE gives higher signal to background ratio as compared to electron beam induced X-ray spectroscopy. As compared to energy dispersive X-ray spectroscopy (EDX), PIXE is 100 times more sensitive. Further PIXE has higher resolution as compared to Rutherford backscattering technique (RBS) [4, 5]. For PIXE analysis the accelerated protons at energy of a few MeV are used to excite the atoms in a sample, which on de-excitation emit characteristics X-rays.

In the PIXE technique the energy of the characteristic X-rays identifies the element, and the number of emitted X-rays with characteristic energy peak is a measure of elemental concentration for specific element [3]. In this technique, accelerated ions beam like alpha particles, protons and other heavy ions have many

applications for the analysis in different research fields like archeological, environmental, biological, geological etc [6-13]. In PIXE technique, the accelerated proton beam is used due to its low bremsstrahlung radiations, high fluorescence yield and X-rays production cross-section [14]. The standardization of PIXE is a process in which materials of already known composition and concentration (Standard Reference Materials (SRM's)) are analyzed for finding the elemental concentration and then compared with the SRM data. For this process it is necessary to set up a reliable calibration of analytical system [15]. The trace elements are those whose concentration are in the parts per billion (ppb) or tenth of parts per million (ppm) in a sample. For standardization of trace elements, accurate knowledge of beam parameters (energy, current and collected charge) and calibration parameters (fluorescence yield, sensitivity, ionization cross section, background radiation and absorption attenuation etc) are required [16].

Thick samples (Titanium based alloy, Copper based alloy and Iron based alloy) were used for the analysis of trace elements in the sample. These samples were bombarded by 2.57 MeV energy proton beam and PIXE spectrum was obtained.

2. EXPERIMENTAL SETUP

Based on our previous work [1] The experimental set up used for PIXE technique is

sketched in the Figure 1. Unlike other spectroscopic techniques, PIXE is based on the physics of characteristic X-rays emitted, not its chemistry [4]. For inducing these X-rays, a high energy proton beam (2.57MeV), accelerated by 5UDH-II tandem Pelletron accelerator excites the atoms of target materials. The excited atoms emit characteristic X-rays in the form of K, L and M X-rays spectrum on its de-excitation [3]. The emitted X-rays are passed through a thin window (Mylar foil 100 μ m made from ultra thin polymer) and detected by a Silicon Drift Detector (SDD) placed at an angle of 45⁰ to the target for good resolution (about 138eV at 5.9KeV energy). The detected signal is amplified by an amplifier. Faraday cup is placed to collect the beam charge and a current integrator records the total collected beam charge during a run. The whole experiment was done under vacuum chamber at a pressure of about 10⁻⁶ torr to reduce the background radiation. After the PIXE spectra is obtained, GUPIXWIN software is used for quantitative analysis of the data [17].

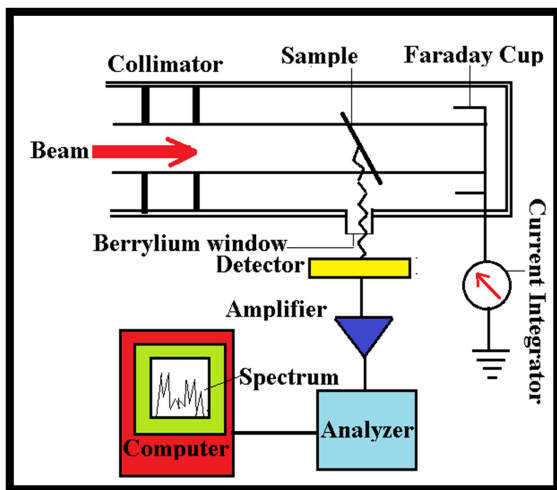


Figure 1: PIXE experimental set up.

The use of GUPIXWIN software is preferred as compared to other software programs due to its good status of databases (e.g. cross-sections, fluorescence yield, Coster–Kronig probabilities, stopping powers and attenuation coefficients etc.) and more elements fitting [18].

3. SAMPLES AND BEAM CHARGE FOR QUANTITATIVE ANALYSIS

For quantitative analysis, PIXE spectrum of each SRM is obtained on different dates for three different values of incident beam charge. Before running the GUPIXWIN code with a proper calibration, we must know elements fit samples

that were used as standard samples were provided by National Institute of Science and Technology (NIST). For two samples NIST SRM 654b (Al 6-V 4) and copper (100% or greater than 99% pure) nine data spectra were acquired by using beam of charge (0.5 μ C, 1 μ C and 1.5 μ C) at three different dates. For NIST SRM 1155 (Cr 18-Ni 12-Mo 2) sample nine data spectra were acquired on a single date to create the parameter files (PAR-files) for the same value of beam charge.

4. EXPERIMENTAL RESULTS

After obtaining the spectrum files, the next important step to get required accuracy is the quantitative analysis by using GUPIXWIN software. Using GUPIXWIN software with suitable parameters, we can determine both elemental composition and concentration present in the target. In GUPIXWIN software “Fixed Matrix solution” method is taken to determine the trace elements in the certified sample. In our experiment, the incident energy beam 2.57 MeV (2570 KeV) and current beam 2nA are used. Position of detector, target and the angle made by detector with target must be fixed during experiment, we only need to determine the instrumental constant (H-value) against the energy data file, i.e. HED file. The determination of instrumental constant (H-value) is the critical point of this research work.

5. INSTRUMENTAL CONSTANT (H-VALUE) Vs ENERGY DATA (HED) FILE

The H quantity is a standardization factor that depends on elements, which translates the peak area of elements of the collected charge into concentration of elements. Initially, in an ideal case, H is taken equal to the solid angle (Ω), that made by the detector in the geometry of the chamber. For our experimental system, H-value is equal to 0.00153 (using formula $\Omega = H = \frac{A}{d^2}$). Practically it slightly depends on the elements atomic number (Z), which produces inaccuracies in detector description and also in data bases [19]. Table 1(a): Combined K HED file

K X-Rays File	
Energy(KeV)	H-Value
1.487	2.5643
1.74	1.95
2.014	0.0145
2.013	0.0077
4.51	0.001147

4.952	0.000975
5.415	0.00096
6.403	0.00098
6.93	0.00415
7.478	0.00097
8.047	0.001288
8.905	0.001118
10.54	0.0049

Table 1(b): L combined HED file

L X-Rays File	
Energy(KeV)	H-Value
1.8066	0.53
2.2932	0.00099
2.3948	1.99
3.444	0.00542
8.3976	0.00053
10.551	0.195

The standard H-value can be determined, either by constant method or by energy dependence method. If geometrical setup is completely accurate then constant method is used and if there is any error then energy dependence method is used. In the later method, HED preparation data file is required. The HED file can be create in note-pad or word-pad with a desired name and need to be saved in the example folder. Further, in the HED file a separate file is present for each (K, L and M) X-ray measurements. It is required to write H-value in front of the high energy (having 100% intensity in the KeV) X-ray values followed by a tab in separate lines for each file (K, L or M X-ray). We need to change the H-values (either increase or decrease) in these files according to the output data result. We have created HED file in two steps.

First, an individual HED file is created for each analyzed sample. Secondly these HED files values are combined to form a combined HED file, which is used for the analysis of all samples. The combined HED files are given in the Table 1(a) and 1(b). These HED files are used to create the PAR-file for required samples. The final concentration obtained is compared with SRM's concentration. Both of concentrations were in good agreement with each other.

6. DATA OF ANALYZED SAMPLES

All samples are calibrated with copper sample using calibration method as previously discussed. The discussion of each sample in detail is given in the following sections.

6.1. NIST SRM 654b (Al 6-V 4) SAMPLE

PIXE curve for NIST SRM 654b (Al 6-V 4) sample is shown in the Figure 2. This curve show both the trace and matrix elements composition and concentration present in sample. As shown in the figure, the identified elements in the sample are Ti, V, Al, Fe, Si, Ni, Cr, Cu and Sn. A parameter (PAR) file of this sample is created, which consist of all information about it. When the PAR-file is completed with a proper sequence, run the gupix


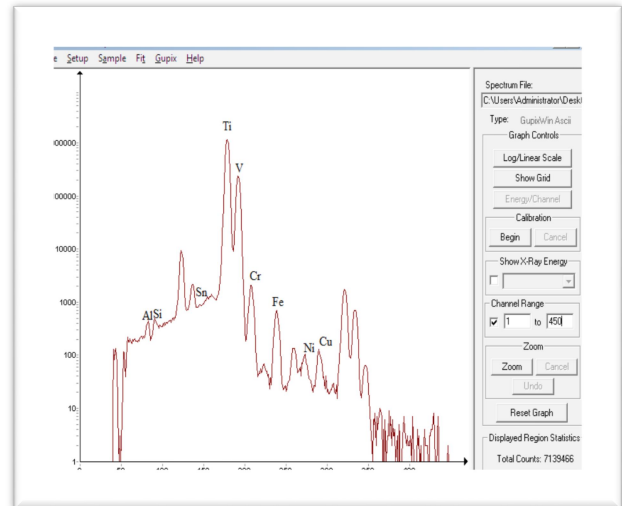


Figure 2: PIXE spectrum for TiV engine to get output data. For TiV sample the concentration of analyzed data, standard data and their standard deviation (σ) are given in Table 2. The dispersion plot (σ) between standard and analyzed data is shown in Figure 3. The standard deviation (σ) values are small, which represent good agreement between measured and already

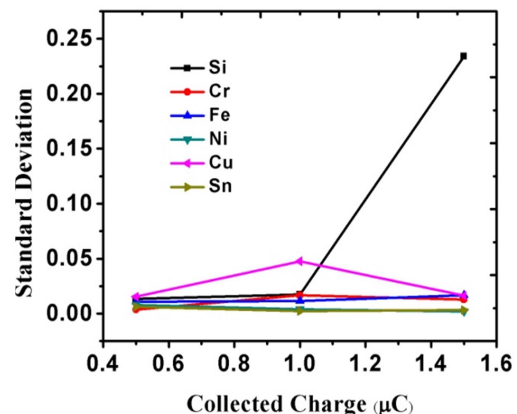

known data.

Figure 3: Dispersion plot for TiV sample.

Table 2: Comparison of measured and actual composition of TiV sample for different beam of charges at different dates.

Charge (μC)	Element	Actual value (a)	Measured value (b_1)	Measured value (b_2)	Measured value (b_3)	Average measured value $\langle b \rangle$	a- $\langle b \rangle$	Standard Deviation (σ)
1.5	Si	0.045%	0.0465%	0.5503%	0.0484%	0.213%	-0.168%	$\pm 0.234\%$
	Cr	0.025%	0.1508%	0.1493%	0.1246%	0.1415%	-0.1165%	$\pm 0.0128\%$
	Fe	0.23%	0.2818%	0.2551%	0.2416%	0.2595%	-0.0295	$\pm 0.0167\%$
	Ni	0.028%	0.0298%	0.027%	0.0251%	0.0273%	0.0007%	$\pm 0.0019\%$
	Cu	0.004%	0.0086%	0.0482%	0.0345%	0.0304%	-0.0264%	$\pm 0.0165\%$
	Sn	0.023%	0.0230%	0.0258%	0.0187%	0.0225%	0.0005%	$\pm 0.0032\%$
1	Si	0.045%	0.0474%	0.0808%	0.0420%	0.0569%	0.0119%	$\pm 0.0174\%$
	Cr	0.025%	0.1850%	0.1563%	0.1454%	0.1621%	-0.1371%	$\pm 0.0166\%$
	Fe	0.23%	0.2753%	0.2666%	0.2512%	0.2643%	-0.0343	$\pm 0.0115\%$
	Ni	0.028%	0.0326%	0.0284%	0.0233%	0.0281%	-0.0001%	$\pm 0.004\%$
	Cu	0.004%	0.0130%	0.0424%	0.0347%	0.03%	-0.026%	$\pm 0.0477\%$
	Sn	0.023%	0.0202%	0.0263%	0.0227%	0.0231%	-0.0001%	$\pm 0.00246\%$
0.5	Si	0.045%	0.0555%	0.0726%	0.0399%	0.056%	-0.011%	$\pm 0.0133\%$
	Cr	0.025%	0.1764%	0.1793%	0.1790%	0.1782%	-0.153%	$\pm 0.0036\%$
	Fe	0.23%	0.2616%	0.2416%	0.2634%	0.2555%	-0.0255	$\pm 0.0107\%$
	Ni	0.028%	0.0386%	0.0277%	0.0195%	0.0286%	-0.0006%	$\pm 0.0078\%$
	Cu	0.004%	0.0091%	0.0426%	0.0412%	0.03096%	-0.027%	$\pm 0.0154\%$
	Sn	0.023%	0.0179%	0.0291%	0.0148%	0.0206%	0.0024%	$\pm 0.0061\%$

6.2. NIST SRM 1155 (Cr 18-Ni 12-Mo 2) Sample.

The PIXE curve for iron base alloy is shown below in the Figure 4.

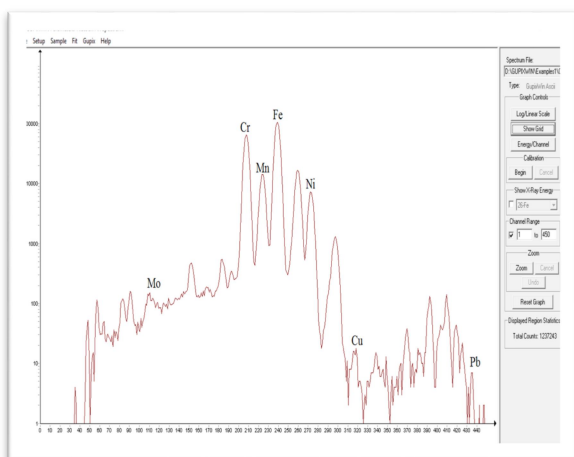


Figure 4: PIXE spectrum for iron base alloy

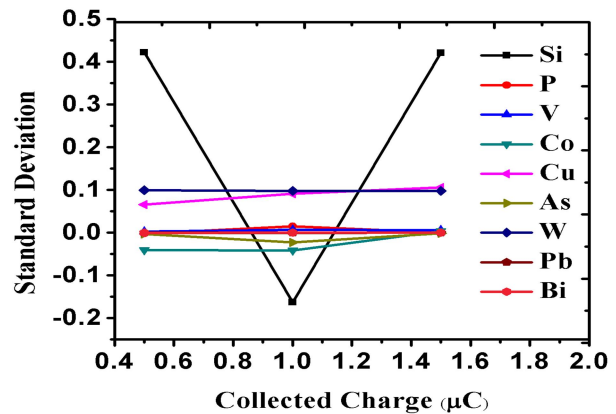


Figure 5: Dispersion plot of iron base alloy

The matrix and trace elements which were identified in this sample are Co, P, Si, V, Cr, Mn, Fe, Ni, Cu, As, Bi and Pb. A parameter (PAR) file is also created using GUPIXWIN software. After running the GUPIXWIN engine, the final output data is recorded as shown in the Table 3. The data was measured on a single date. Difference between measured data and already known data were shown in Table 3. The dispersion plot for this sample is shown in the Figure 5.

The difference between concentration of measured and already known data values are negligibly small.

Table 1: Comparison of measured and actual composition of iron sample for different beam of charge

Charge (μC)	Element	Actual value (a)	Measured value (b)	(a-b)
1.5	Si	0.5093%	0.0889%	0.4204%
	P	0.0200%	0.0201%	-0.0001%
	V	0.050%	0.0437%	0.0063%
	Co	0.109%	0.1055 %	0.0035%
	Cu	0.175%	0.0690%	0.106%
	As	0.01067%	0.0107%	-0.00173%
	W	0.11%	0.0124%	0.0976%
	Pb	0.001 %	0.0011%	-0.0001%
	Bi	0.0004%	0.0002%	0.0002%
1	Si	0.5093%	0.6723%	-0.163%
	P	0.0200%	0.0047%	0.0153%
	V	0.050%	0.0430%	0.007%
	Co	0.109%	0.5256%	-0.0417%
	Cu	0.175%	0.0840%	0.091%
	As	0.01067%	0.0335%	-0.0228%
	W	0.11%	0.0124%	0.0976%
	Pb	0.001 %	0.0012%	-0.0002%
	Bi	0.0004%	0.0002 %	0.0002%
0.5	Si	0.5093%	0.0802%	0.4219%
	P	0.0200%	0.0224%	-0.0024%
	V	0.050%	0.0470%	0.003%
	Co	0.109%	0.1497%	-0.0407%
	Cu	0.175%	0.1093%	0.0657%
	As	0.01067%	0.0137%	-0.0030%
	W	0.11%	0.0108%	0.0992%
	Pb	0.001 %	0.0011%	-0.0001%
	Bi	0.0004%	0.0006%	-0.0002%

6.3 COPPER (Cu) SAMPLE

By using GUPIXWIN software, PIXE curve obtained for copper sample is shown in the Figure 6. The matrix and trace elements identified in this sample

are Cu, Ti and Sr. The major element is Cu, with

concentration of 99% or greater, while the other two (Ti and Sr) are trace elements (concentration less than 1%).

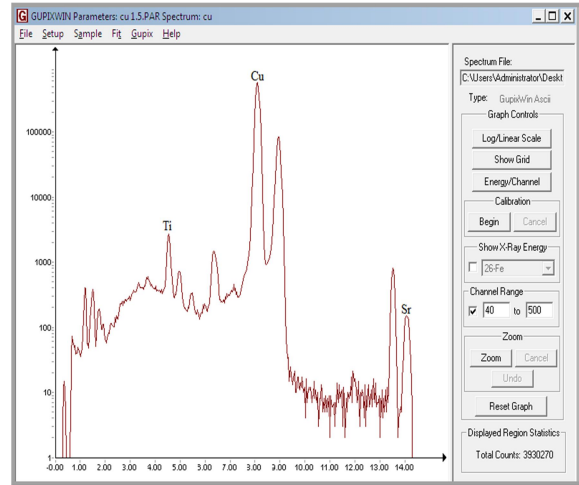


Figure 6:PIXE spectrum for Copper Sample.

Standard data for Ti and Sr were unknown. The concentrations found for trace elements (Ti and Sr in Cu sample for different beam charges at different dates and have small standard deviation A PAR-file is created for Cu sample which include all information about it.

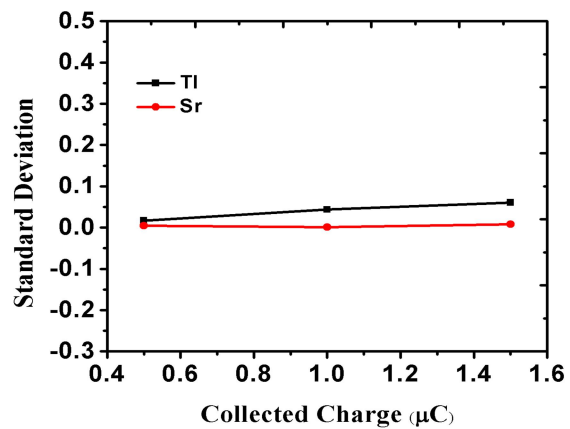


Figure 7: Dispersion plot of Cu sample

The Table 4 show the concentration of trace element identified in copper sample (more than 99% concentration) at different dates and the results obtained are found to be in close agreement. The standard deviation for Cu sample is given in Table 4 and plotted in plotted in Figure 7. It is clear from the Figure that the difference of measured and already known data is negligibly small.

Table 2: Comparison of measured and actual composition of copper sample for different beam of charges at different date

Charge (μC)	Element	Actual Value (a)	Measured Value (b ₁)	Measured Value (b ₂)	Measured Value (b ₃)	Average 	Standard Deviation (σ)
1.5	Ti	-----	0.22%	0.08%	0.08%	0.12%	$\pm 0.0602\%$
	Sr	-----	0.07%	0.08%	0.0614%	0.07%	$\pm 0.0079\%$
1	Ti	-----	0.16%	0.07%	0.07%	0.10%	$\pm 0.0435\%$
	Sr	-----	0.0795%	0.0989%	0.0754%	0.08%	$\pm 0.0012\%$
0.5	Ti	-----	0.11%	0.07%	0.0832%	0.08%	$\pm 0.0168\%$
	Sr	-----	0.0860%	0.0914%	0.09%	0.0905%	$\pm 0.0048\%$

7. CONCLUSIONS

In our previous paper [1], we have shown the calibration of PIXE technique for thick samples. So the main objective of this research work is to collaborate the PIXE technique for trace element analysis. Three samples (two NIST and one local) were analyzed using PIXE technique. Trace elements composition of measured data and standard data of known composition were in close agreement. In TiV sample the maximum standard deviation (σ) for Si is 0.0234%, Cr 0.0166%, Fe 0.0167%, Ni 0.0078%, Cu 0.047% and for Sn is 0.0061%. Similarly, in Cu sample, the maximum deviation for Ti is 0.0602% and for Sr is 0.0079%. The standard deviation values for all trace elements in TiV and Cu were very small (less than 1%). In iron base alloy sample, the maximum difference is very small (less than 1%) for each trace element. The very small (less than 1%) differences in elemental composition of measured and known values show that PIXE is versatile, feasible and reliable technique for trace elements analysis.

8. BIBLIOGRAPHY

- [1] Shad Ali, Johar Zeb, Abdul Ahad, Ishfaq Ahmad, M. Haneef, Jehan Akbar, "Standardization of proton-induced x-ray emission technique for analysis of thick samples," Chinese Physics B, vol. 24, no. 090601, 2015, DOI:10.1088/1674-1056/24/9/090601
- [2] Johansson S A E and Johansson T B 1976 Nuclear Instrument and Methods **137** 473
- [3] Ene A, Popescu I V, Stihl C, Gheboianu A, Pantelica A and Petre C 2009 Rom. Jour. Phy **55** 806
- [4] Lin E K, Wang C W, Yu Y C, Liu T Y, Tau T P, Chiou J W 1997 Chinese Journal of Physics **35** 6
- [5] Johansson S A E, Campbell J L, and Malmqvist K G 1995 Particle-Induced X-ray Emission Spectroscopy (New York: John Wiley & Sons) pp.147-412
- [6] Kabir M H 2007 Particle Induced X-ray Emission (PIXE) Setup and Quantitative Elemental Analysis (Ph.D. Dissertation) (Kochi: University of Technology Japan)
- [7] Vijayan, Behera V, Ramamurthy S N, Puri V S, Sanjiv, Shahi, Singh J S and Nirmal 1997 Nuclear Instrument and Methods **26** 65
- [8] Bandhu H K et al 2000 Nuclear Instrument and Methods B **160** 126
- [9] Ene V, Popescu I V and Stihl C 2010 International Journal of Physics **55** 806
- [10] Kennedy V J, Aughtly A, Variar K M, Magudapathy P, Panchapakesan, Vijayan V and Chakravorty 2002 **254** 351
- [11] Nayak P K, Das D, Chintalapudi S N, Singh P, Acharya S, Vijayan V and Chakravorty V J 2002 Radioanalytical Nuclear Chemistry **254** 351
- [12] Nayak P K, Das D, Vijayan V, Singh P and Chakravorty V 2001 Nuclear Instrument and Methods B **184** 649
- [13] Hajivaliei M et al 1999 Nuclear Instrument and Methods **150** 645
- [14] Rath A K, Acharya G S and Vijayan V 1993 Journal Methods of Physics **18** 1
- [15] Chakravorty A, Selvaraj S, Sudarshan M, Dutta R K, Ghugre S S and Chintalapudi S 2000 Nuclear Instrument and Methods B **170** 156
- [16] Kasahara M, Ogiwara G, Yamamoto K, Park J H and Takahashi 1994 Int. Journal of PIXE **4** 155
- [17] Cohen D D, Siegle R, Orlic I, and Stelcer E 2001 Long term Accuracy and precision of PIXE and PIGE measurements for thin and thick sample analysis, June 8-12, 2001, Menai, Australia, p. 1

- [18] Govil M 2001 Current Science **80** 1542
- [19] Khak P M, Kouhpar M K, Hajivaliei and Khademi F 2013 Mediterranean Archaeology and Archaeometry **13** 1
- [20] Shariff A, Bulow K, Elfman M, Kristiansson P, Malmqvist K and Pallon J 2002 Nuclear Instruments and Methods in Physics Research B **189** 131