

Commissioning of low particle flux for proton beams at MedAustron

Felix Ulrich-Pur^{a,*}, Laurids Adler^c, Thomas Bergauer^a, Alexander Burkert^b, Andrea De Franco^c, Greta Guidoboni^c, Albert Hirtl^b, Christian Irmler^a, Stefanie Kaser^a, Sebastian Nowak^c, Florian Pitters^a, Mauro Pivi^c, Dale Prokopovich^c, Claus Schmitzer^c, Alexander Wastl^c

^aAustrian Academy of Sciences, Institute of High Energy Physics (HEPHY), Nikolsdorfer Gasse 18, 1050 Wien, Austria

^bTU Wien, Atominstitut, Stadionallee 2, 1020 Wien, Austria

^cMedAustron, Marie-Curie-Straße 5, 2700, Wiener Neustadt, Austria

Abstract

MedAustron is a synchrotron-based particle therapy centre located in Wiener Neustadt, Austria. It features three irradiation rooms for particle therapy, where proton beams with energies up to 252.7 MeV and carbon ions of up to 402.8 MeV/u are available for cancer treatment. In addition to the treatment rooms, MedAustron features a unique beamline exclusively for non-clinical research (NCR). This research beamline is also commissioned for proton energies up to 800 MeV, while available carbon ion energies correspond to the ones available in the clinical treatment rooms.

Based on the requirements for particle therapy, all irradiation rooms offer particle rates of up to 10^9 particles/s for protons and 10^7 particles/s for carbon ions. However, for research purposes, lower particle fluxes are required and were therefore commissioned for the NCR beamline. Three particle flux settings with particle rates ranging from $\approx 2.4 \times 10^3$ particles/s to $\approx 5.2 \times 10^6$ particles/s were established for seven proton energies below 252.7 MeV. In addition to the particle rate, the spot sizes and beam energies were measured for these settings. Furthermore, three low flux settings for 800 MeV protons with particle rates ranging from $\approx 2 \times 10^3$ particles/s to $\approx 1.3 \times 10^6$ particles/s were commissioned. Since the commissioned low flux settings are in a regime well below the limits of the available standard beam diagnostics, setting up the beam under these new operational conditions entirely relied on the use of external detectors. Furthermore, a beam position measurement based alignment without using the standard beam profile monitors was performed for 800 MeV protons.

Keywords: MedAustron, low particle flux, intensity reduction, synchrotron

1. Introduction

MedAustron is a particle therapy and research centre located in Wiener Neustadt, Austria. It features four irradiation rooms (IR1-IR4) with one beamline (IR1) exclusively dedicated to research [1]. Protons up to 800 MeV and carbon ions of up to 402.8 MeV/u can be delivered to this beamline. The design of the MedAustron accelerator complex was optimized by medical requirements [2], resulting in available particle rates in the order of $\approx 10^8 - 10^9$ particles/s for medical use. However, experiments in non-clinical research often require measuring the interaction of single particles [3], which leads to completely different beam and detector requirements compared to those for particle therapy. Based on the accelerator design [4], the maximum beam intensity could be theoretically lowered by a factor of 650, leading to a minimum particle rate of $\approx 10^6$ particles/s. But due to the rising demand for lower particle rates far below the design specifications [4] ($\approx 10^3 - 10^6$ particles/s), low flux settings for proton beams were commissioned in the research beamline at MedAustron. However, the beam for the required low flux settings had to be set up “blindly”, as the conventional beam diagnostic elements were designed for use

with therapeutic beam intensities down to $\approx 10^7$ particles/s. Therefore, the lower limit of the available clinical intensities ($\approx 10^8$ particles/s), where beam instrumentation still could be used, was taken as the initial set-up for the low flux settings. From this point, the initial set-up was then extrapolated down to the low flux levels without the use of the standard beam instrumentation. Instead, to set up and monitor the low flux beam, external detectors were placed at the isocentre of IR1.

2. Materials and methods

2.1. MedAustron accelerator

Figure 1 shows the MedAustron layout. As primary ions, H_3^+ or $^{12}C^{4+}$ are generated in two electron cyclotron resonance (ECR) ion sources at 8 keV/u [5] located in the low energy beam transfer line (LEBT). The ions are then transported to the radio frequency quadrupole (RFQ) followed by an interdigital H-mode linear accelerator (LINAC), which accelerates the particles up to 7 MeV/u. After the LINAC, a carbon stripping foil is used to strip off remaining electrons and convert the ions to H^+ or C^{6+} . Those ions are then transported to the synchrotron by the medium energy beam transfer line (MEBT) and injected via a multi-turn injection scheme. After injection, the beam is bunched and accelerated by a radio frequency cavity. Inside the

*Corresponding author

Email address: felix.ulrich-pur@oeaw.ac.at (Felix Ulrich-Pur)

synchrotron ring (77 m circumference), protons are accelerated up to 800 MeV and carbon ions up to 402.8 MeV/u. Via a third order slow resonant extraction [6] the accelerated beam is then extracted towards one of the four irradiation rooms through the High Energy Beam Transfer line (HEBT).

2.2. Particle flux reduction methods

The flexible design of the MedAustron accelerator complex allows to expand the possible range of particle rates. Different methods were combined in order to reduce the particle rate per area. These are:

1. Reducing the number of particles injected into the synchrotron using the electrostatic fast deflector (EFE) or a degrader.
2. Extending the extraction time into the HEBT to reduce the number of particles per second by changing the betatron core ramping.
3. Scraping the beam into the chopper dump or increasing the transverse beam spot size in the HEBT so it exceeds the size of the vacuum tube. The increased size of the beam leads to a loss of particles in the HEBT and to a larger spot size at the isocentre, resulting in a reduced number of particles per unit area.

It is important to notice that not all low flux settings mentioned above were applied for proton beams in the medical energy range (between 62.4 and 252.7 MeV) and at 800 MeV. The increased transverse beam size method was not used in the case of 800 MeV to avoid environmental irradiation at that energy. In the following sections, a detailed description of the accelerator elements and the method used for flux reduction is given.

2.2.1. Electrostatic fast deflector

The EFE consists of two electrodes installed in the LEBT line. If a high voltage is applied (5.5 kV [7]), the beam is deflected onto a cylindrical faraday cup, otherwise the beam can pass through towards the LINAC. The main purpose of the EFE is actually to adjust the pulse length: a typical injection scheme foresees 30-50 μ s pulses for the phase space painting in the synchrotron, and it can be triggered with a 1 μ s resolution, therefore providing a well controllable measure to reduce the beam intensity arriving in the main ring.

2.2.2. Degrader

In the MEBT section, a pepper pot like device (degrader) is used to adapt the number of particles injected into the ring [8]. The beam intensity from the LINAC can be reduced to 10 %, 20 % or 50 % of the nominal value.

2.2.3. Betatron core

The slow extraction mechanism used at MedAustron is the betatron core third order driven resonant extraction [9]. The use of the betatron core features the advantage of keeping all lattice elements constant during extraction and therefore the beam properties in the irradiation room.

At first, the beam is set slightly below the synchrotron's design energy as defined by the dipole's magnetic field. This causes the bunch to circulate along the dispersion orbit (the synchrotron RF and the resonant sextupole are located in a dispersion free region). During extraction, the machine is tuned such that the extraction occurs when the zero-betatron amplitude particle has reached the designed energy and thus circulates in the centre of the beam pipe.

In preparation of the extraction, an RF jump is performed in order to create a more homogeneous $\Delta p/p$ distribution. Then the RF is turned off and a coasting beam is obtained after a few ms. The resonant sextupole is ramped while the beam is still positioned safely away from the resonant tune. Finally, the beam is slowly driven into the resonance by inductively accelerating it with the betatron core. The acceleration effect is given by the flux change $\Delta\Phi = C \cdot B\rho \cdot \Delta p/p$ [10], where C is the machine circumference and $B\rho$ is the magnetic rigidity of the particle. As the beam is accelerated into the resonance, the large emittance particles become unstable first: their amplitude grows until they reach the electrostatic septum that deflects them into the extraction line. The betatron core ramping can be adjusted to extend the extraction time and reduce the particle rate.

2.2.4. Chopper

The chopper is a system located in the HEBT line (see Figure 1) comprised of four dipoles powered in series that allow a fast (less than 300 μ s) switching of the beam (controlled switch on/off of the beam during normal operations and for emergency interruptions) [5]. Only when the dipoles are powered, the beam can circumvent the dump and be transmitted along the extraction line (see Figure 2). For low flux beam generation, the chopper has been used to control the transmission through the dump. The dipole strength was adjusted such that part of the beam was scraped on the dump and only a reduced portion of the extracted beam from the synchrotron reached the room. This method has been applied to protons at 800 MeV only.

2.2.5. Quadrupole magnets to IR1

The last 7 quadrupoles to IR1 - the irradiation room for non-clinical research purposes - were used to blow up the transverse (to the direction of motion) size of the beam and therefore reduce the particle flux at the monitor location. The modification of quadrupole settings implies a so-called optics change of the transfer line and the method will be referenced as optics adjustment in the following. Figure 3 shows a comparison of the β_x and β_y functions along the transfer line to IR1. The β_x and β_y functions represent the transverse beam envelopes which are related to the beam size in the transverse planes. For the same beam, a larger beta function implies a larger beam size. The quadrupole setting (Optics5) is used in case of the nominal setting in contrast to new settings with increased transverse beam size, Optics3 and Optics4, respectively.

2.3. Low flux settings

As mentioned in section 2.2, different methods for implementing low flux were applied depending on the beam energy.

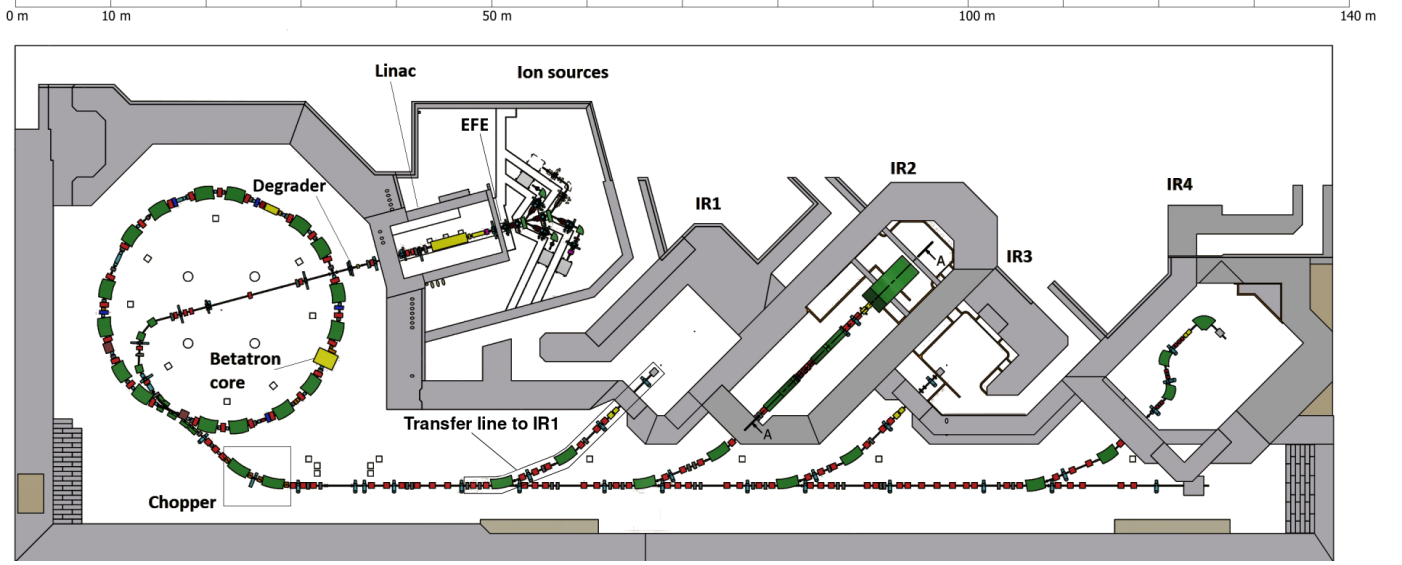


Figure 1: Layout of the MedAustron accelerator complex where the key elements for low flux beam generation are shown. The (small) chopper magnets are located around the two HEBT dipoles in the box on the left. In the transfer line to IR1, the last seven quadrupoles are used to increase the transverse size of the beam for flux reduction.

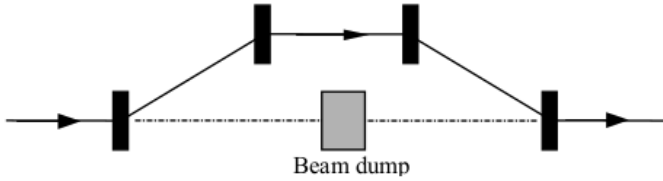


Figure 2: Schematic of the chopper system. The black rectangles represent the four dipoles and the solid arrows show the beam path during normal operation. When the chopper is off, the beam hits the dump (grey box).

Table 1 summarizes the combination of accelerator elements to achieve different particle rates (particles/s), called “Set. MED” in case of medical energy range (between 62.4 MeV and 252.7 MeV) or “Set. p800” for protons at 800 MeV. For each energy range, three settings (“high”, “medium”, “low”) were commissioned. The order of magnitude of the obtained mean particle rate per setting defines the name of each setting and will be given in units of Hz. Setting “high”, “medium” and “low” correspond to a mean particle rate of $O(\text{MHz})$, $O(100 \text{ kHz})$ and $O(\text{kHz})$, respectively.

2.4. Detectors used for beam characterization

The experimental set-ups for measuring the rate, spot sizes and beam energies are described in detail in the following section.

2.4.1. Rate monitor

Since the current beam diagnostics of MedAustron were not designed for low flux beams, a dedicated rate monitor (RaMon) to achieve single particle counting at low particle rates

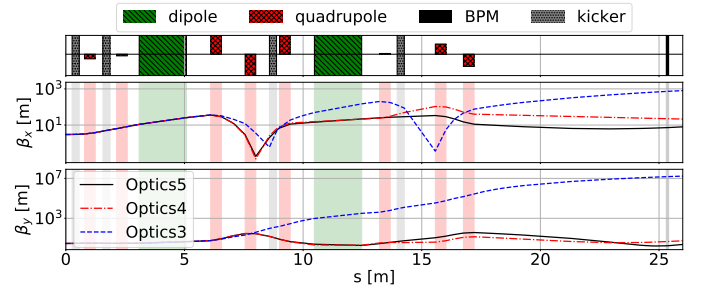


Figure 3: On the upper side of the figure, details of the transfer line are shown. Quadrupole magnets are highlighted in red, dipole and kicker magnets in green and grey, respectively, and beam position monitors (BPM) in black. The β functions along the transfer line to IR1 are displayed on the central (horizontal, x-plane) and bottom plots (vertical, y-plane) on a logarithmic scale. The nominal optics is shown in black, being very similar to Optics4 displayed in red. Optics3 is in blue and much larger in both planes compared to the other settings.

(< 50 MHz) was built for the rate measurements.

For this purpose, two EJ228 plastic scintillators from Eljen [11] were chosen as particle counters, since they have a rise and decay time of a few ns. Each of the scintillators has a total volume of $50 \times 50 \times 10 \text{ mm}^3$. As light guide, PMMA fish tail light-guides ($50 \times 10 \text{ mm}^2$) were connected to the photocathode ($\varnothing = 8 \text{ mm}$) of a Hamamatsu H10721-210 photosensor [12]. Optical grease and double-sided optical duct tape were used for the assembly of each scintillator. The Hamamatsu photomultipliers (PMTs) were powered and readout by the AIDA 2020 trigger and logic unit (TLU) [13]. A coincident signal of both PMTs defines a single particle count. Each coincident event was assigned with a timestamp using the internal 160 MHz clock of the TLU. Every 100 ms, the number of coincidences and re-

Set. MED	high	medium	low
Appr. rate	4.7 MHz	350 kHz	3 kHz
EFE	1 μ s	1 μ s	1 μ s
Degrader	10 %	10 %	10 %
Betatron	Nom.	10 % Nom.	10 % Nom.
Optics	Optics4	Optics5	Optics3
H. size σ_x	3.9-8.6 mm	3.2-8.5 mm	> 5 cm
V. size σ_y	2.8-8.8 mm	2.4-7.6 mm	> 5 cm
Set. p800	high	medium	low
Appr. rate	1.25 MHz	200 kHz	2 kHz
EFE	30 μ s	10 μ s	30 μ s
Degrader	20 %	20 %	20 %
Betatron	Nom.	Nom.	12 % Nom.
Chopper	60.5 % Nom.	60.5 % Nom.	60.5 % Nom.
H. size σ_x ¹	8.2 mm	9.4 mm	8 mm
Set. NOM	nominal		
Appr. rate	200 MHz		
EFE	30 μ s		
Degrader	10 %		
Betatron	Nom.		
Optics	Optics5		
H. size σ_x	2.8-8.7 mm		
V. size σ_y	2.9-8.9 mm		

¹ Only the horizontal spot size σ_x was measured for 800 MeV (see section 3.2).

Table 1: The commissioned low flux settings are listed. For proton beams in the medical energy range, the settings include the combination of EFE, degrader and/or betatron core and optics adjustments. In case of protons at 800 MeV, the chopper is used instead of optics adjustment. The approximate mean rate and the horizontal and vertical spot size (σ_x , σ_y) are given for each setting. For comparison, the used nominal setting is also shown.

spective timestamps were processed by the Data Acquisition System (DAQ). For this purpose, the EUDAQ2 framework [14] was chosen since an interface for the AIDA TLU is already implemented in this software framework. An online beam monitor was developed using EUDAQ2 and QT [15] to calculate and display the particle rates. 30-50 spills were used per setting to calculate the average rate. A schematic overview of the RaMon system can be seen in Figure 4 on the right.

2.4.2. Spot size measurement

Depending on the particle flux, two different detectors for measuring the spot size at the IR1 isocentre were used. For all low flux settings, a double sided silicon strip detector (DSSD), which has already been used for single particle tracking at MedAustron [3], was used. The n-substrate based DSSD is 300 μ m thick and has an active area of 2.56×5.12 cm². Each side features 512 orthogonal AC coupled strips with a pitch of 100 μ m on the n-side (x -coordinate) and 50 μ m on the p-side (y -coordinate). The sensor itself is readout via 4 APV25 chips [16] on each side (128 strips/APV, Figure 4). A custom readout system [3], which was originally developed for the Belle II Silicon Vertex Detector [17], was used. The readout was trig-

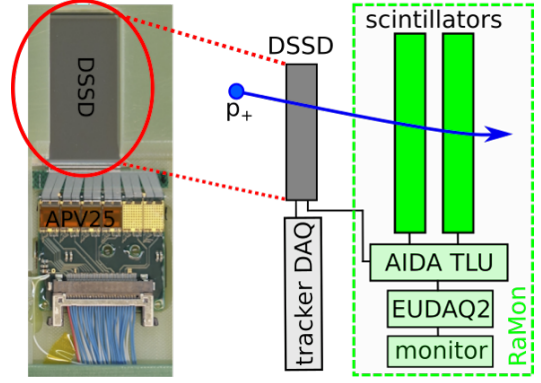


Figure 4: Schematic overview of the spot size measurement for low flux protons using the DSSD and the RaMon system. On the left, an image of the p-side (y -coordinate) of the DSSD is shown. The DSSD strips are connected to 4 APV25 chips mounted on hybrid boards which are readout with a custom DAQ system. The RaMon plastic scintillators are connected to the AIDA2020 TLU and are used to trigger single particle events.

gered using the RaMon plastic scintillators, which were placed behind the DSSD. A schematic overview of the experimental setup is shown in Figure 4.

Due to the limited particle rate capability of the RaMon system (< 50 MHz), the Lynx[®] detector (IBA-Dosimetry, Schwarzenbruck) was used for the nominal setting (≈ 200 MHz). This detector has already been used for spot size and position verification at clinical beam intensities during commissioning in the other treatment rooms [18], [19]. The Lynx[®] detector itself is a gadolinium-based plastic scintillating screen with an active area of 30×30 cm² and a thickness of 0.4 mm. Via a mirror, the scintillating screen is coupled to a CCD camera with 600×600 pixel and a pixel size of 4.25×4.25 μ m² [20]. All components of the Lynx[®] detector are contained in a light-tight box with $360 \times 370 \times 600$ mm³. A schematic overview of the Lynx[®] is depicted in Figure 5.

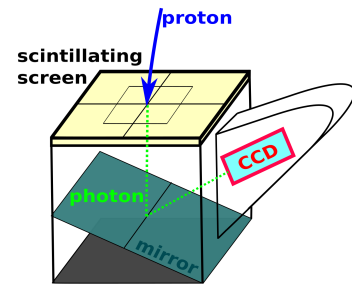


Figure 5: Schematic overview of the spot size measurement with the Lynx[®] detector. The interaction position of an incident particle on the 30×30 cm² scintillating screen is mapped onto the 600×600 pixel of a CCD camera via a mirror.

Whenever a particle deposits energy in the scintillating screen, the deposited energy is converted into green light (540 nm), which is then reflected towards a CCD camera via a mirror. The amount of light reaching the CCD camera can be modified by an aperture collimator system (iris) to avoid satura-

tion. Furthermore, the exposure time for a single image can be adapted (≤ 180 s). According to [20], the effective spatial resolution of the Lynx[®] detector is 0.5 mm. Even though the Lynx[®] detector could also have been used for particle rates down to $O(100 \text{ kHz})$, the DSSD system was used instead for all low flux settings due to its superior spatial resolution.

2.4.3. Energy measurement

In addition to the spot size and rate, the influence of each low flux setting on the beam energy was studied by measuring the range in water for different beam energies and flux settings. The range R in a particular material can be obtained using the continuous-slowing-down approximation (CSDA). In order to calculate the CSDA range for a proton with initial energy E_0 , the reciprocal total stopping power $S(E)$ has to be integrated between the proton's initial state $E = E_0$ and stationary state $E=0$.

$$R_{\text{CSDA}} = \int_{E_0}^0 \left(-\frac{dE}{dx} \right)^{-1} dE = \int_{E_0}^0 \frac{dE}{S(E)} \quad (1)$$

For a homogeneous material, simplified, parameterized versions of equation (1) ([21], [22]) can be used to convert the measured range R_{CSDA} into the initial beam energy E_0 . Therefore, the range obtained in a residual range detector is usually expressed in water equivalent thickness (WET[23]). If the traversed material in the range detector consists of n different, homogeneous materials, with total stopping powers $S_i(E)$ ($i \in [1, n]$) and the traversed pathlength per material is t_i , equation (1) can be rewritten as follows

$$R_{\text{CSDA}} = \underbrace{\int_{E_0}^{E_1} \frac{dE}{S_1(E)}}_{\hat{=t_1}} + \underbrace{\int_{E_1}^{E_2} \frac{dE}{S_2(E)}}_{\hat{=t_2}} + \dots + \underbrace{\int_{E_{n-1}}^0 \frac{dE}{S_n(E)}}_{\hat{=t_n}} \quad (2)$$

$$\Rightarrow R_{\text{CSDA}} = \sum_i^n t_i. \quad (3)$$

In order to convert the measured CSDA range into range in water, the water equivalent thickness WET_i of each traversed pathlength t_i leading to R_{CSDA} has to be calculated according to [23]

$$\text{WET}_i \approx t_i \frac{S_i(E)}{S_{\text{H}_2\text{O}}(E)} \Rightarrow R_{\text{CSDA, WET}} = \sum_i^n \text{WET}_i. \quad (4)$$

For the following measurements, the total stopping power in water $S_{\text{H}_2\text{O}}(E)$, as well as the total stopping power for each traversed material $S_i(E)$ were taken from the NIST database [24] to obtain the final residual range in water $R_{\text{CSDA, WET}}$.

For particle rates above $O(100 \text{ kHz})$, the residual range in water was measured directly using the PTW PEAKFINDER[™] (PTW, Freiburg, Germany). The PTW PEAKFINDER[™] (Figure 6) is a height-adjustable water column with a diameter of 114 mm and two ionization chambers to measure the depth dose profile of a particle beam in water. The water column itself consists of two interconnected water-filled bellows (absorber and reservoir). The PTW TM34082 is used as a reference chamber (RC) in front of the absorber water column and the PTW TM34080 as a field chamber (FC) between the absorber and reference water column. The position of the FC as well as the length of each

bellow is adjusted via a servo motor (stepsize $\geq 10 \mu\text{m}$). To sample the depth-dose profile, the absorbed dose to water [25] in the FC is measured relative to the absorbed dose to water in the RC at different positions inside the water column (340 mm movable range).

Close to the maximum (Bragg peak), the depth-dose profile was sampled in 0.3 mm steps. The range in water was then measured at the 20 % distal fall-off relative to the Bragg peak for each energy and compared to the ranges from the NIST PSTAR database.

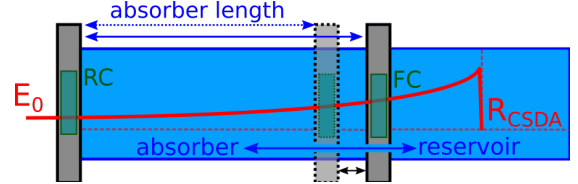


Figure 6: The PTW PEAKFINDER[™] consists of two interconnected, water-filled bellows (absorber and reservoir) and 2 ionization chambers (RC, FC). The depth-dose profile in water was obtained by measuring the dose to water at different positions of the FC in the water column. From the sampled depth-dose profiles, the range in water was derived for different flux settings and beam energies.

Due to a low signal-to-noise ratio in the ionization chambers below particle rates of $O(100 \text{ kHz})$, no depth-dose profiles could be obtained with the PTW PEAKFINDER[™]. Instead, a range telescope, formerly developed by the TERA foundation [26], was used to measure the particle range in water. Fast silicon photomultipliers (Hamamatsu MPPC S10362-11-050C [27]) coupled to 38 plastic scintillator slices with a water equivalent thickness of 3.6 mm and an active area of $30 \times 30 \text{ cm}^2$ each allow range measurements of single particles. The readout of a single particle event is triggered by the coincident signal of the Ra-Mon plastic scintillators placed in front of the range telescope. The measured range in the telescope is then converted to residual range in water according to equation (4). As described in [3] and [28], the range telescope suffered from severe voltage instabilities, therefore the mainboard and the readout software were completely redesigned and replaced prior to the measurements.

2.5. Additional challenges

For normal (medical) operational conditions, the beam position and profile in the HEBT can be measured with scintillating fibre hodoscopes (SFX) [29]. Due to the low beam current of the low flux proton beams, the SFX detectors in the extraction line do not produce enough signal to measure a beam profile, effectively rendering the transfer line “blind” to the beam position until the isocentre. Therefore, under these new operational conditions, the verification of the beam position and beam angle at the isocentre completely relied on external detectors. In order to measure the beam angle at the isocentre, three additional DSSDs were placed behind the first DSSD at the isocentre at equal distances (10 cm). The measured angles were found to be negligible for all settings.

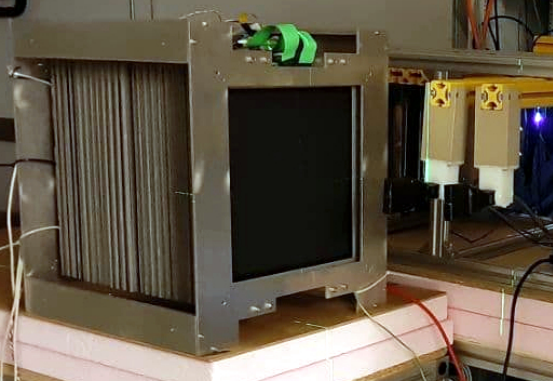


Figure 7: Range measurement for low flux protons at MedAustron using the TERA range telescope with two plastic scintillators in front as triggers.

2.5.1. Steering of the 800 MeV low flux beam

During commissioning of the 800 MeV low flux beam, apart from a reduced strength of the chopper dipoles, the setpoints of all other elements in the extraction line were initially left unchanged compared to the nominal flux beam. During measurement of the spot sizes with the DSSD, an offset of the beam in the isocentre of 12.5 mm horizontally and -13 mm vertically could be observed. In order to correct the beam position at the isocentre while not changing the beam angle, a second position measurement upstream would be required. Since that proved to be impossible due to the low beam current, an adapted steering strategy was employed, where the corrector kick response matrix of the beam position [30], using the last two horizontal and vertical correctors upstream the IR1 isocenter, was measured for the nominal flux beam on the SFX close to the last quadrupole of the beamline and at the isocentre. The kick response matrix describes the change in beam position as a function of a difference in applied kick angle $\Delta\Theta_{x,j}$ by a corrector magnet. The elements of the kick response matrix are defined as (exemplary for the horizontal plane)[30]

$$C_{ij}^{xx} = \frac{\Delta x_i}{\Delta\Theta_{x,j}} [\text{m/rad}], \quad (5)$$

where Δx_i is the change of the horizontal beam position at the beam profile monitor BPM_i [m] and $\Delta\Theta_{x,j}$ is the change of the horizontal kick angle of corrector j [rad]. After obtaining the kick response matrix via measurement of the position response, the required change in corrector kick angles $\Delta\vec{\Theta}_{x,\text{needed}}$ to correct the beam position on all beam profile monitors to its target position (reference trajectory) can be calculated by inverting the kick response matrix

$$\begin{pmatrix} \Delta\vec{\Theta}_{x,\text{needed}} \\ \Delta\vec{\Theta}_{y,\text{needed}} \end{pmatrix} = \begin{pmatrix} C^{xx} & C^{xy} \\ C^{yx} & C^{yy} \end{pmatrix}^{-1} \cdot \begin{pmatrix} \Delta\vec{x}_{\text{target}} \\ \Delta\vec{y}_{\text{target}} \end{pmatrix}. \quad (6)$$

Where $\Delta\vec{x}_{\text{target}}$ is the vector of the difference between the (horizontal) beam position and the target beam position on all beam profile monitors. The inverted kick response matrix of the nominal flux beam was then applied to the measured offset of the

low flux beam, resulting in new corrector strengths compensating the offset while still minimizing the beam angle at the isocentre. Without being able to measure the low flux beam with the BPM_i in the HEBT, this beam position measurement based alignment was performed using the DSSDs at the isocentre only. With this method, the 800 MeV low flux beam was aligned prior to the rate and spot size measurements with 800 MeV protons.

3. Results

As described in Section 2.3, three different flux settings could be commissioned for beam energies below 252.7 MeV and three different settings for 800 MeV. For details of the settings, the reader is referred to Table 1.

3.1. Rate measurement

The particle rates for each setting were measured with the RaMon system for seven different beam energies below 252.7 MeV and three different settings for 800 MeV. The combined results are shown in Figure 8.

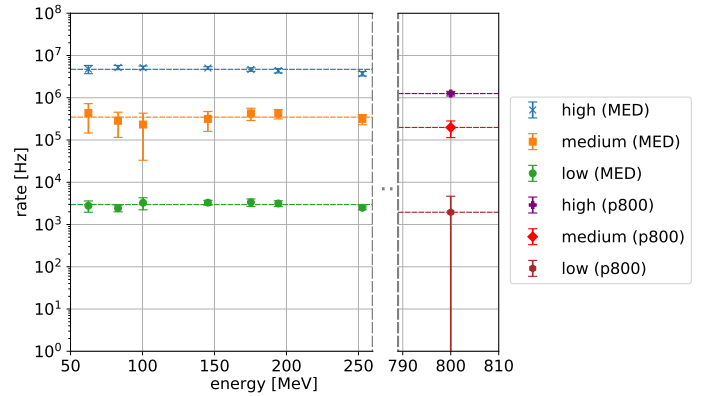


Figure 8: Resulting low flux particle rates measured with the RaMon system. Three different flux settings per energy could be commissioned. The settings for 800 MeV protons differed from the medical settings. The dashed horizontal lines represent the mean particle rate per setting, which also defines the name of the settings.

The mean rate for the lowest medical flux setting, called “low (MED)”, ranges from 2.43 – 3.35 kHz. Depending on the particle energy, a mean rate between 232 – 435 kHz was obtained for setting “medium (MED)”. For setting “high (MED)”, mean particle rates between 3.73 – 5.21 MHz were measured. The measured rates for the 800 MeV settings showed similar results compared to the medical settings. A mean particle rate of 1.95 kHz for setting “low (p800)”, 198 kHz for setting “medium (p800)” and 1.25 MHz for setting “high (p800)” were obtained. Setting “low (p800)” showed more significant rate fluctuations per spill compared to the other 800 MeV low flux settings. Figure 9 depicts the spill structure for setting “low (p800)” in comparison to setting “high (p800)”.

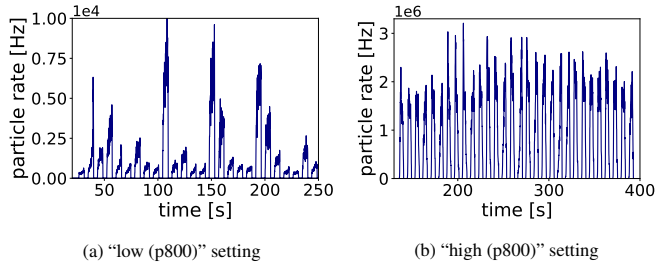


Figure 9: Comparison of the spill structure of low flux protons for 800 MeV protons. Setting “low (p800)” shows larger fluctuations in measured particles per spill compared to setting “high (p800)”.

3.2. Spot size measurement

The Lynx[®] detector was used to measure the spot sizes for the nominal setting (≈ 200 MHz) in a prior measurement [31]. A 2D Gaussian was fitted to the obtained 2D intensity profile distribution. The resulting full widths at half maximum (FWHM) were compared to spot sizes of the low flux medical settings, which were measured with the DSSD. The beam spots for the “high (MED)” and “medium (MED)” settings also showed a Gaussian-shaped distribution. The measured horizontal spot size is depicted in Figure 10 and the vertical spot sizes in Figure 11.

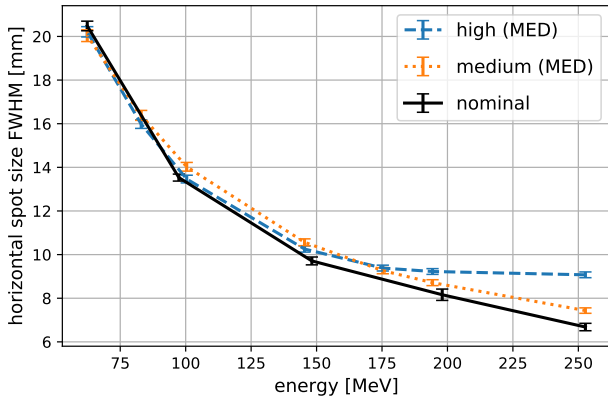


Figure 10: Horizontal (x-axis) beam spot size for different proton energies and particle fluxes.

Similar vertical spot sizes compared to the nominal flux setting were obtained for the “high (MED)” setting. Except for the horizontal spot size for “high (MED)” at 252.7 MeV, the absolute difference in FWHM compared to the nominal flux was always less than ≈ 1 mm FWHM. The vertical spot size of the “medium (MED)” setting was slightly smaller compared to the nominal flux ($\approx 1 - 2$ mm Δ FWHM). On the other hand, the absolute difference in FWHM of the “medium (MED)” setting compared to the nominal flux did not exceed ≈ 0.5 mm FWHM.

The beam profiles of the “low (MED)” setting were also measured with the DSSD. Figure 12 shows the measured vertical beam profile of the “low (MED)” setting for 252.7 MeV protons as an example.

It is apparent that the shape of the beam is always bigger than

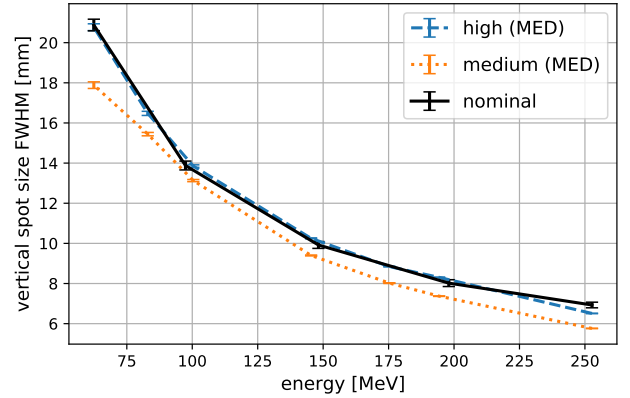


Figure 11: Vertical (y-axis) beam spot size for different proton energies and particle fluxes.

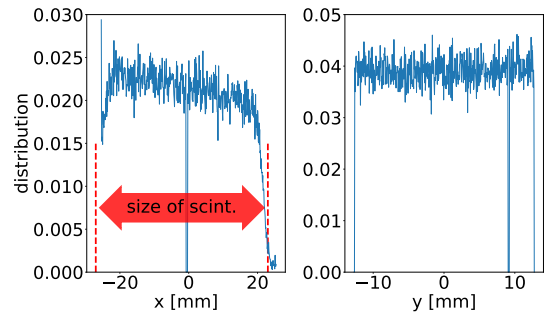


Figure 12: Horizontal (left) and vertical (right) beam profile of 252.7 MeV protons for setting “low (MED)” measured with a DSSD. The spot size cannot be determined since the beam is bigger than the trigger scintillators used for the DSSD. The beam profile of the horizontal beam close to the trigger edges is truncated due to the size of the trigger scintillators.

the 5×5 cm² plastic scintillators. Only events passing both scintillators are able to trigger an event in the DSSD, which leads to a truncated beam profile at the edges of the plastic scintillators ($\approx \pm 25$ mm). This also means that the measured rate for this setting can only be related to the area of the plastic scintillators, which has to be taken into account when using this setting. The broadening of the beam is a direct result of the use of the magnetic quadrupoles in the HEBT. Due to the geometry of the trigger scintillators and the DSSD, only part of the beam profile distribution could be measured. However, the beam distribution appears to be relatively uniform inside the measured area (5×2.56 cm²). This was also observed for the horizontal beam profile.

The spot size for 800 MeV protons was also measured with the DSSD. Prior to the spot size measurement, the horizontal and vertical offset of the 800 MeV beam was corrected as described in Section 2.5. In contrast to the medical settings, the spot sizes for this energy were not measured in the isocentre, since they were recorded simultaneously to the rate measurement with the RaMon in the isocentre. For this purpose, the DSSD was placed 16 cm downstream the RaMon. Thus, the measured spot size is increased due to multiple Coulomb scattering in the plastic scintillators of the RaMon setup in front of

the DSSD. The beam profiles for all three 800 MeV low flux proton settings are depicted in Figure 13.

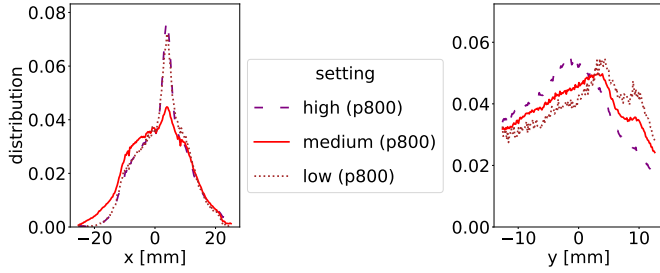


Figure 13: Horizontal (left) and vertical (right) beam profile of 800 MeV protons for setting “low (p800)”, “medium (p800)” and “high (p800)” measured with a DSSD. The beam profiles do not follow a Gaussian distribution. The vertical beam is slightly larger than the vertical side of the DSSD.

From Figure 13, it is apparent that the beam profiles for all 800 MeV settings do not follow a Gaussian distribution. A second peak near the centre of the spot could be detected for all settings. The vertical beam profile for all settings is slightly larger than the p-side (vertical axis) of the DSSD. On the other hand, the horizontal beam profiles are smaller than the n-side (horizontal axis) of the DSSD and appear centred. In order to estimate the spot size, the sample standard deviation σ_x of the horizontal beam profile was calculated. It resulted in 8.2 mm for the “high (p800)”, 9.4 mm for the “medium (p800)” and 8 mm for the “low (p800)” setting.

3.3. Energy measurement

The effect of the particle flux reduction method on the mean energy was studied by measuring the range in water for different flux settings and different beam energies. The measured range in water was then compared to the NIST PSTAR database. Because of an insufficiently low signal-to-noise ratio in the ionization chambers of the PTW PEAKFINDER™, the TERA range telescope was used to measure the range in water for the “low (MED)” setting at 83 and 100.4 MeV. For all other flux settings, including the nominal flux, the PTW PEAKFINDER™ could be used. As can be seen in Figure 14, no significant change in the measured range was observed for any flux setting. The measured ranges are also in good agreement with the ranges obtained from the NIST database.

In contrast to the rate and spot size measurements, the energy measurements were recorded prior to a timing optimization of the accelerator cycle to reduce dead time. Even though no change in the particle energy was expected, the measurement for the “medium (MED)” setting at 252.7 MeV was repeated after the timing optimization. As expected, the same range was obtained.

4. Summary and Outlook

In total, six low flux settings for protons were commissioned at MedAustron and can now be used in the non-clinical research room. Because the particle rate of the commissioned

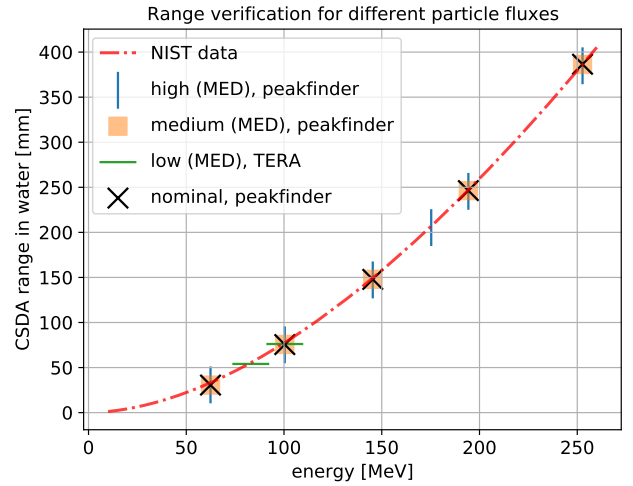


Figure 14: Measured range in water for different proton energies and flux settings. The obtained ranges were compared to data from the NIST PSTAR database.

settings was well below the design specifications of the accelerator and beam diagnostic elements, the beam for the required settings had to be set up “blindly” without the standard beam instrumentation, only using external detectors. The beam energy, particle rate and spot size were measured for three settings at seven different beam energies ranging from 62.4 to 252.7 MeV. Furthermore, three low flux settings were commissioned for 800 MeV.

A dedicated rate monitor system (RaMon), consisting of two plastic scintillators, was developed (Figure 4) for the particle rate measurements. For each energy, a stable setting with a mean particle rate of \mathcal{O} (kHz), another with \mathcal{O} (100 kHz) and a setting with \mathcal{O} (MHz) was commissioned.

The beam profiles were measured using the Lynx® for higher particle rates in the clinical energy range and a DSSD for all low flux settings. Setting “high (MED)” and “medium (MED)” showed similar beam profiles compared to the nominal (≈ 200 MHz) flux setting. However, a slightly smaller vertical spot size was obtained for the setting “medium (MED)” ($\approx 1 - 2$ mm Δ FWHM)). A difference in the width of the measured horizontal beam profiles larger than 1 mm FWHM was only observed for the “high (MED)” setting at 252.7 MeV. The effect of the magnetic quadrupoles used for the “low (MED)” setting resulted in a large beam spread, hence the full beam profile could not be measured for this setting. The measured rate for setting “low (MED)” can only be defined on the area as large as the plastic scintillators (50×50 mm²). The initially measured beam profiles for all 800 MeV settings showed a vertical and horizontal offset in the isocentre. Therefore, a beam position measurement based alignment had to be performed prior to the spot size and rate measurements for all p800 settings. This was done solely with external detectors (DSSDs) since the standard beam diagnostic elements in the HEBT were not designed for particle rates well below the therapeutic intensities ($\approx 10^3 - 10^6$ particles/s). For all three 800 MeV settings, the

measured beam profiles did not follow a Gaussian distribution. To estimate the spot size, the sample standard deviation was calculated for the horizontal beam profile. Depending on the setting, a sample standard deviation ranging from 8 – 9.4 mm was obtained. In order to measure the full vertical beam profile for 800 MeV, bigger detectors, able to handle low particle fluxes, are required.

The range in water for different beam energies and flux settings was measured using the PTW PEAKFINDER™ and a range telescope based on plastic scintillator slices coupled to SiPMs. A comparison of the measured ranges in water for different particle energies and particle fluxes yielded that the particle flux reduction methods had no significant impact on the particle energy.

Acknowledgements

This project has received funding from the European Union's Horizon 2020 research and innovation programme under the Marie Skłodowska-Curie grant agreement No 675265 and from the Austrian Research Promotion Agency (FFG), grant number 869878.

References

- [1] M. Benedikt (Ed.), Physics opportunities at MedAustron : white book; 1. ed., Schriftenreihe der Technischen Universität Wien, Univ. of Technology, Vienna, 2009, Literaturverz. S. 83 – 100. URL <https://repository.gsi.de/record/67546>
- [2] M. Benedikt, J. Gutleber, M. Palm, et al., Overview of the MedAustron design and technology choices (2010) 3 p. URL <https://cds.cern.ch/record/1272158>
- [3] F. Ulrich-Pur, T. Bergauer, A. Burkert, et al., Imaging with protons at MedAustron, Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment 978 (2020) 164407. doi:10.1016/j.nima.2020.164407.
- [4] L. Badano, M. Benedikt, P. Bryant, et al., Proton-ion medical machine study (PIMMS) part II, CERN/PS 2000-07(DR) (2000).
- [5] M. Benedikt, A. Wrulich, MedAustron—Project overview and status, The European Physical Journal Plus 126 (7) (2011) 69. doi:10.1140/epjp/i2011-11069-9.
- [6] M. Benedikt, P. Bryant, M. Pullia, A new concept for the control of a slow-extracted beam in a line with rotational optics: Part II, Nuclear Instruments and Methods in Physics Research Section A: Accelerators, Spectrometers, Detectors and Associated Equipment 430 (2) (1999) 523 – 533. doi:10.1016/S0168-9002(99)00207-7.
- [7] J. Borburgh, T. Fowler, T. Kramer, et al., Design of Electrostatic Septa and Fast Deflector for MedAustron (CERN-ATS-2011-270) (2011) 4 p. URL <http://cds.cern.ch/record/1407937>
- [8] L. Adler, A. De Franco, F. Farinon, et al., Synchrotron Emittance Analysis Procedure at MedAustron, Proceedings of the 7th Int. Beam Instrumentation Conf. IBIC2018 (2019) China. doi:10.18429/JACoW-IBIC2018-WEPC02.
- [9] M. G. Pullia, et al., Betatron Core Driven Slow Extraction at CNAO and MedAustron, in: Proc. of International Particle Accelerator Conference (IPAC'16), Busan, Korea, May 8-13, 2016, no. 7 in International Particle Accelerator Conference, JACoW, Geneva, Switzerland, 2016, pp. 1330–1333. doi:10.18429/JACoW-IPAC2016-TUPMR037.
- [10] L. Badano, M. Benedikt, P. Bryant, et al., Proton-ion medical machine study (PIMMS) part I, CERN/PS 99-010(DI) (1999) 105.
- [11] ELJEN Technology, FAST TIMING EJ-228, EJ-230 (2021).
- [12] Hamamatsu, Hamamatsu photosensor module H10721-210 (2016).
- [13] P. Baesso, D. Cussans, J. Goldstein, The AIDA-2020 TLU: a flexible trigger logic unit for test beam facilities, Journal of Instrumentation 14 (09) (2019) P09019–P09019. doi:10.1088/1748-0221/14/09/p09019.
- [14] Y. Liu, M. Amjad, P. Baesso, et al., EUDAQ2—A flexible data acquisition software framework for common test beams, Journal of Instrumentation 14 (10) (2019) P10033–P10033. doi:10.1088/1748-0221/14/10/p10033.
- [15] The QT framework, <https://www.qt.io/>, accessed: 2021-04-20.
- [16] M. French, et al., Design and results from the APV25, a deep sub-micron CMOS front-end chip for the CMS tracker, Nucl. Instr. and Meth. A 466 (2) (2001) 359 – 365. doi:10.1016/S0168-9002(01)00589-7.
- [17] R. Thalmeier, et al., The Belle II SVD data readout system, Nucl. Instr. and Meth. A 845 (2017) 633 – 638. doi:10.1016/j.nima.2016.05.104.
- [18] M. Stock, D. Georg, A. Ableitinger, et al., The technological basis for adaptive ion beam therapy at MedAustron: Status and outlook, Zeitschrift für Medizinische Physik 28 (3) (2018) 196–210. doi:<https://doi.org/10.1016/j.zemedi.2017.09.007>.
- [19] L. Grevillot, J. O. Moreno, V. Letellier, et al., Clinical implementation and commissioning of the MedAustron particle therapy accelerator for non-isocentric scanned proton beam treatments, Medical Physics 47 (2) (2019) 380–392. doi:10.1002/mp.13928.
- [20] S. Russo, A. Mirandola, S. Molinelli, et al., Characterization of a commercial scintillation detector for 2-D dosimetry in scanned proton and carbon ion beams, Physica Medica 34 (2017) 48 – 54. doi:10.1016/j.ejmp.2017.01.011.
- [21] W. Bragg, R. Kleeman, On the alpha particles of radium, and their loss of range in passing through various atoms and molecules, The London, Edinburgh, and Dublin Philosophical Magazine and Journal of Science 10 (57) (1905) 318–340. doi:10.1080/14786440509463378.
- [22] W. Donahue, W. D. Newhauser, J. F. Ziegler, Analytical model for ion stopping power and range in the therapeutic energy interval for beams of hydrogen and heavier ions, Physics in Medicine and Biology 61 (17) (2016) 6570–6584. doi:10.1088/0031-9155/61/17/6570.
- [23] R. Zhang, W. D. Newhauser, Calculation of water equivalent thickness of materials of arbitrary density, elemental composition and thickness in proton beam irradiation, Physics in Medicine and Biology 54 (6) (2009) 1383–1395. doi:10.1088/0031-9155/54/6/001.
- [24] M. Berger, J. Coursey, M. Zucker, et al., Stopping-power and range tables for electrons, protons, and helium ions, NIST Standard Reference Database 124 (2017). doi:10.18434/T4NC7P.
- [25] Absorbed Dose Determination in External Beam Radiotherapy, no. 398 in Technical Reports Series, INTERNATIONAL ATOMIC ENERGY AGENCY, Vienna, 2001.
- [26] M. Bucciantonio, et al., Development of a fast proton range radiography system for quality assurance in hadrontherapy, Nucl. Instr. and Meth. A 732 (2013) 564 – 567. doi:10.1016/j.nima.2013.05.110.
- [27] Hamamatsu, MPPC S10362-11 Datasheet (2013).
- [28] M. Bucciantonio, Development of an advanced Proton Range Radiography system for hadrontherapy, Dissertation, Philosophisch-naturwissenschaftliche Fakultät der Universität Bern (2015). URL https://swisscovery.sls.ch/permalink/41SLSP_NETWORK/1ufb5t2/alma991084316869705501
- [29] M. Feurstein, A. Gyorgy, A. Kerschbaum, et al., Overview of the Beam Diagnostics in the MedAustron Accelerator: Design Choices and Test Beam Commissioning, Conf. Proc. C 1205201 (2012) 774–776.
- [30] H. Wiedemann, Particle Accelerator Physics, Springer International Publishing, 2015. doi:10.1007/978-3-319-18317-6.
- [31] H. Fuchs, B. Knäusel, P. Kuess, et al., Abnahmebericht nichtklinischer Bestrahlungsraum IR1 (unpublished), Tech. rep., MedAustron (2016).